




SENIOR DISSECTING MICROSCOPE MD-2

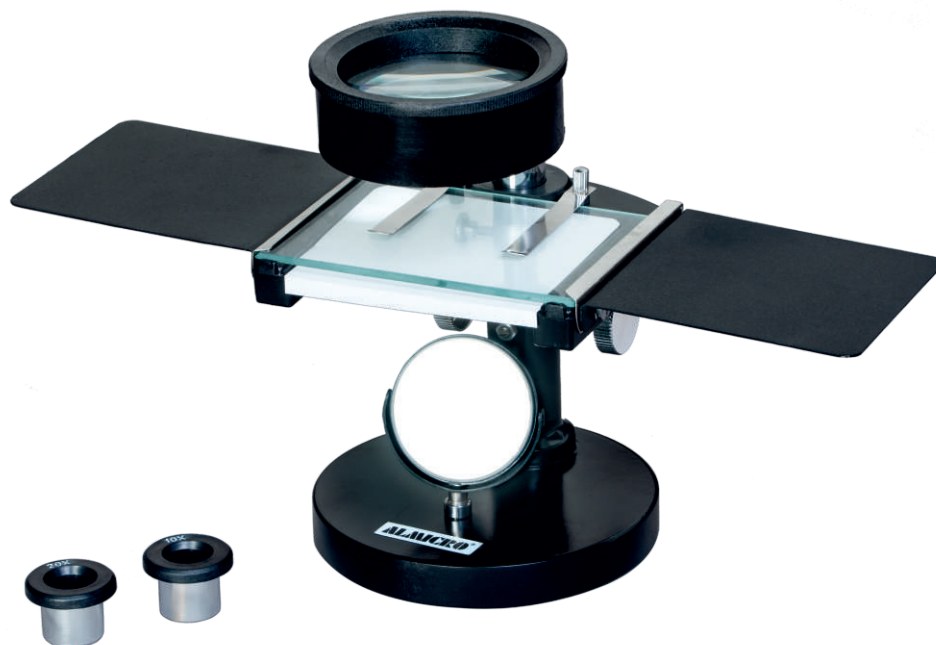

MD-2 Senior Dissecting Microscope has heavy round base with a precisely designed body. Sensitive focussing is done by rack-pinion movement. A joint revolving arm is provided for moving magnifying lens overfull area of the stage. It is easily replaceable by a special bulls eye lens (83mm dia.) to increase the field of vision.

A piano-concave mirror reflects sufficient day light and artificial light on the object. Detachable hand rests on both sides of stage provided for a study support during dissection. Provided with 10x & 20x eye pieces and special bulls eye lens (83mm dia.).

Packed in wooden carrying cabinet.

Packaging: Inner Thermocol/Styrofoam Packing and Outer Cardboard Box.